

Liquid handling

PUMPS AND ULTRA-HIGH PRESSURE VALVES

Diluter/dispensers, M series

- Full liquid handling functionality
- Self-priming
- No syringes
- Largest volume range available
- Easy to use Wizard format does away with math problems and charts



Cheminert M Series diluter/dispensers simplify the sample preparation process for dispensing and diluting liquids. A user-friendly Wizard format eliminates all the math calculations and charts usually associated with diluting and dispensing applications. Just enter the dilution ratio and the final volume, and the correct volume is calculated and automatically dispensed for each ratio.

For multiple dispenses, simply enter the volume and the number of dispense repetitions, and the Wizard calculates the total volume to be aspirated. It's that easy!

The diluter/dispenser is built around a patented syringe-free, bi-directional, positive displacement pump. This design approach gives the largest volume range available, and eliminates the inconvenience of having to change and refill syringes.

M10 Diluter/Dispensers

	Prod No	Price
10 nL - 10 mL - M10 diluter/dispenser	CD10-4841-M1A	\$3500

M50 Diluter/Dispensers

	Prod No	Price
50 µL - 50 mL - M50 diluter/dispenser	CD50-8182-M2A	\$3600

Additional Features

"Smart" hand probe

The hand probe signals the operator when an aspirate or dispense step is completed. The unique design also allows the use of fixed or disposable probe tips, as well as other accessories.

Program memory

Up to 100 programs can be permanently stored.

Multi-solvent option

A multiposition stream selection valve can be easily integrated with the pump for multi-solvent applications.

Printer option

Print out methods, sequential steps, time/date/operator stamp, titration, and tubing volume values.

Applications

- Simple dispensing of reagents using the manual dispense mode
- Micro dispensing in microplates and genomic arrays
- Dilutions for AA, ICP, GC, and HPLC samples
- Serial dilutions for all samples
- Multi-sample and reagent additions, micro-plates, tube to plate, tube to tube
- Small and large volume dispensing of reagents
- Titrations

Liquid handling pumps, M series

CE

The Cheminert® M Series liquid handling pump is a syringe-free pump capable of delivering a bidirectional flow over six orders of magnitude.

The M Series is a positive displacement pump, which means that it is self-priming and tolerant of any gas which may find its way into the fluid lines. There is no separate fill cycle, and the capacity is unlimited.

RS-232 and RS-485 communication protocols are incorporated into the microprocessor-driven controller. (USB interface requires an adapter.) The included software package controls flow rates, flow direction, and metered volumes.



Operating principle

At the core of the pump is a polymeric rotor housing four 1/8" diameter pistons in sapphire cylinders. As the microstepper motor turns the rotor, the pistons float on a stationary cam; at any given moment, one piston is filling, one is dispensing, and the other two are in transit between the fill and dispense positions.

SPECS

	M6	M50
Continuous minimum dispense	100 nL	50 µl
Continuous maximum dispense	5 ml/min	25 ml/min
Maximum back pressure	100 psi	100 psi
Gravimetric precision		
for 125 µl	0.5%	0.8%
for 1.25 ml	0.05%	0.1%
Pump internal volume (µl)	100 ± 2 µl	625 ± 10 µl

Applications

- Flow cytometry, cell and drug perfusion
- HTS and robotic systems
- Infusion and micro-dialysis
- Micro diluters/dispensers for nL to mL range applications
- Micro liquid transfers (nL) for micro arrays
- Microtiter plate dispensing using multiposition valves

TECH TIP

Use a standoff assembly if the motor must be separated from the pump head. Standoffs are available in lengths of 2", 3", 4", and 6".

M6 pumps

10 nL - 10 mL

	Prod No	Price
M6 pump with:		
Controller and stepper motor	CP2-4841-100M1	\$2740
Stepper motor (no controller)	CP2-4841-100SM	2060
M6 pump only	CP2-4841-100D	1540

M50 pumps

50 µl - 50 mL

	Prod No	Price
M50 pump with:		
Controller and stepper motor	CP3-8182-625M2	\$2920
Stepper motor (no controller)	CP3-8182-625SM2	2230
M50 pump only	CP3-8182-625D	1710

Accessories and replacement parts

	Prod No	Price
Pump motor		
M6	CP-DSM	\$525
M50	CP-DSM2	550
Controller, MicroLynx-4	CP-CM1-P	600
Standoff assembly		
2"	2SOAMP	50
3"	3SOAMP	50
4"	4SOAMP	50
6"	6SOAMP	50

Ultra-high pressure injector system

LIQUID HANDLING

40,000 psi ultra-high pressure injector system

The VICI 40K injector is comprised of six miniature air actuated needle valves, plumbed to simulate the flow path of a conventional rotor/stator injector. An integral controller sends the on/off positioning signals to each valve, coordinating them to perform load, inject, and flush functions.


There are three methods for sending positioning commands to the injector:

- Manual control with the pushbuttons on the controller
- Laboratory computer via serial port communication
- Contact closure inputs


For more information, contact our technical department.

Ultra-high pressure injector system

Prod No	Price
SPSS40	\$4000

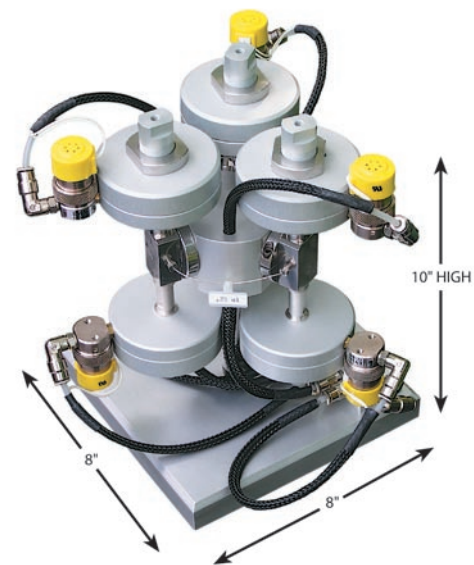


Ultra-high pressure injector system schematic

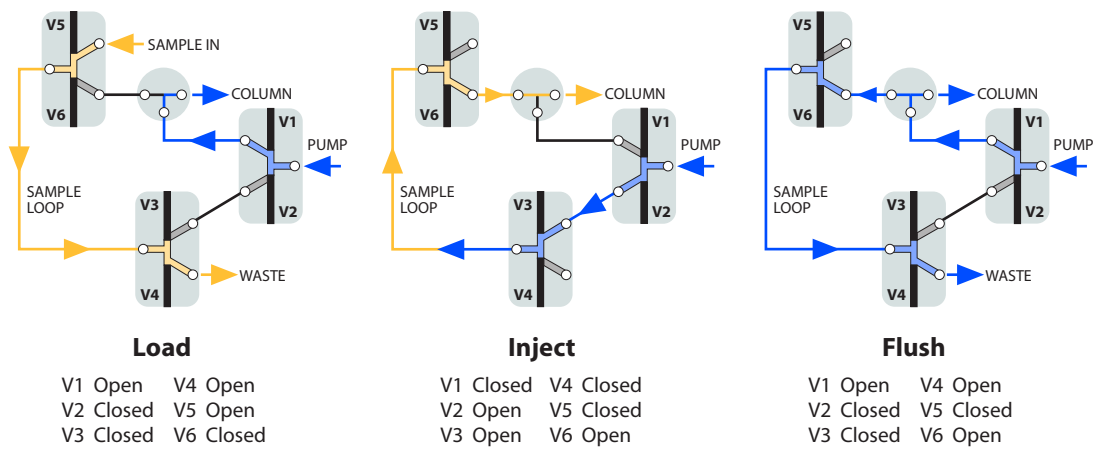

40,000 psi ultra-high pressure valves

The ultra-high pressure valves that are the heart of our SPSS40 (*previous page*) are now available individually, in 1/16", 1/32", and 360 micron versions. There are three types – a two port on/off valve, a dual on/off valve, and a 3-way prime/purge valve. (*See page 213 for flowpath schematics.*) The dual on/off configuration has two individually controlled outlets with a common inlet (or vice versa), emulating a rotary three way valve.


Implementation requires a single three-way solenoid: application of 50 psi opens the valve; venting the air allows the spring to return the valve to the closed position. A fitting for 1/8" air supply tubing is included; two fittings are included for dual valves. (*Fitting: prod no EAOR21, page 219*)

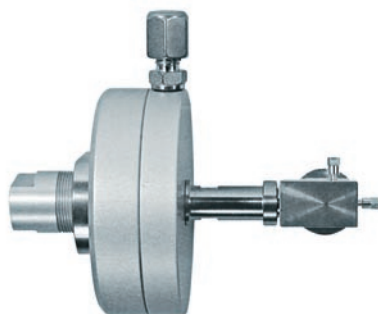

On/off valve
360 µm ZDV fittings

On/off valves

SPECS		Fitting size	Bore	Prod No	Price
Temp	Pressure				
50°C	40,000 psi	360 µm	0.15 mm	ASFVO40K360	\$840
		1/32"	0.15 mm	ASFVO40K.5	760
		1/16"	0.15 mm	ASFVO40K1	715

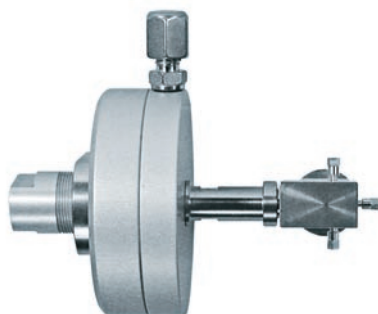


Prime/purge valve
360 µm ZDV fittings

Prime/purge valves


SPECS		Fitting size	Bore	Prod No	Price
Temp	Pressure				
50°C	40,000 psi	360 µm	0.15 mm	ASFV40K360	\$880
		1/32"	0.15 mm	ASFV40K.5	790
		1/16"	0.15 mm	ASFV40K1	745

Dual on/off valves

SPECS		Fitting size	Bore	Prod No	Price
Temp	Pressure				
50°C	40,000 psi	360 µm	0.15 mm	ASFVOD40K360	\$1460
		1/32"	0.15 mm	ASFVOD40K.5	1340
		1/16"	0.15 mm	ASFVOD40K1	1260

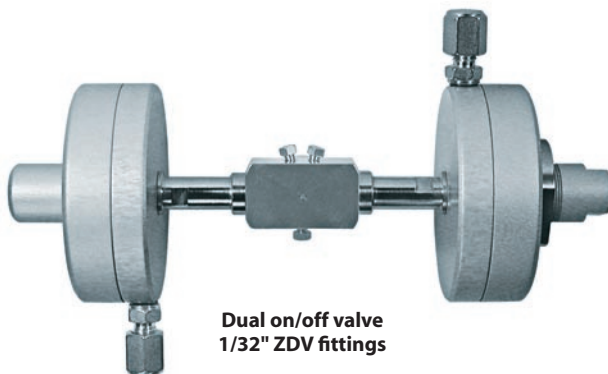

Dual on/off valve
1/32" ZDV fittings

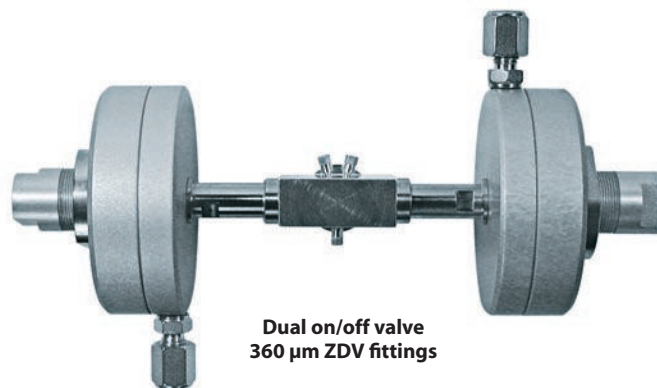

Dual on/off valve
360 µm ZDV fittings



Equipamientos Analíticos S.A

DESDE 1989 | BRINDANDO SOLUCIONES